

PAL-RI

for measuring refractive index

Cat#3850

1.3

134 1 342

425

25

20 c

Dear Distributor,

We are proud to announce that on this date, October 21st, 2005, we will launch our newest digital "Pocket" refractometer PAL-RI.

This long-awaited digital "Pocket" refractometer is a compact sized instrument that measures the refractive indices of liquid samples and it is equipped with all the other great features of our conventional PAL series refractometers except the temperature compensation.

for the chemical and pharmaceutical related industries.



By pressing the START key, after placing some droplets of the sample on the prism, the refractive index and the temperature of the measurement will be displayed continuously like an electric newsboard.

Equipped with digital display, it is lightweight, small in size, and water-resistant too. It also features the E.L.I. technology ensuring that excessive light will not affect the measurement.

Zero setting (calibration) is easy using water only. And a wide measurement range of refractive index of 1.3306 to 1.5284 will make it an ideal model

| Specifications | | | |
|--------------------------------|---|------------------------|---|
| Model | PAL-RI | Cat.No. | 3850 |
| Measurement Range | Refractive index 1.3306 to 1.5284 | Resolution | Refractive index 0.0001 |
| Measurement Accuracy | Refractive index ±0.0003 (water at 20°C) | Ambient Temperature | 10 to 40°C |
| Measurement Temperature | 5 to 45°C (resolution 1°C) | Sample Volume | 0.3 ml |
| Measurement Time | 3 seconds | Power supply | 2 x AAA Batteries |
| International Protection class | IP65 Dust-tight and Protected against water jets. | Dimensions & weight | 55(W) x 31(D) x 109(H)mm 100g (Main unit only) |

For more details and prices, please call your ATAGO contact now!

ATAGO CO.,LTD.

32-10, Honcho, Itabashi-ku, Tokyo, 173-0001, Japan Email: export@atago.net Web: http://www.atago.net/ TEL: 81-3-3964-6156 FAX: 81-3-3964-6137

ATAGO U.S.A.,Inc.

12011 Bel-Red Road, Suite 101, Bellevue, WA 98005, U.S.A. E-mail: customerservice@atago-usa.com TEL: 1-425-637-2107 TOLL FREE: 1-877-ATAGO-USA FAX: 1-425-637-2110

ATAGO INDIA Instrumente Pvt. Ltd.

Office 2 Nilkanth Smruti, shimpoli Rood, Borivali (West), Mumbai-400 092 India

TEL: 91-22-2833-8038 / 8076 FAX: 91-22-2899-8164 E-mail: customerservice@atago-india.com

Options

- · 30% Sucrose (±0.03%): RE-110030
- · 40% Sucrose (±0.04%): RE-110040
- · 50% Sucrose (±0.05%): RE-110050
- · 60% Sucrose (±0.05%): RE-110060
- · Daylight plate for small volume samplings : RE-2391-50M
- * Calibration Certificate : Contact an ATAGO representative for details.

· Daylight plate : RE-2315-61M


Design Registration No.000379326-0001, 000379326-0002 (EU) ZL 2005 3 0116403.4,5 (China), D111526 (Taiwan), D555, 549 (U.S.A.),

1255763, 1255764, 1255765, 1255766, 1255767 (Japan) Patent Granted in countries around the world.

Refractometer MASTER-2M





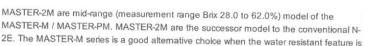

2E. The MASTER-M series is a good alternative choice when the water resistant feature is not a priority and the Temperature Compensation will be carried out manually. Not equipped with the Automatic Temperature Compensation feature. The MASTER-2M is made of metal, and the MASTER-2PM is plastic. The specifications are

exactly the same. For salty and acidic samples, choose the plastic body, MASTER-2PM. *Calibration with the 30% Sucrose Solution (RE-110030), 40% Sucrose Solution (RE-110040), 50% Sucrose Solution (RE-110050) or 60% Sucrose Solution (RE-110060)

60 Brix 30





Specifications

| Model | MASTER-2M |
|-------------------------|----------------------|
| Cat.No. | 2323 |
| Scale | Brix |
| Measurement Range | Brix: 28.0 to 62.0% |
| Minimum Scale | Brix: 0.2% |
| Measurement Accuracy | Brix: ±0.2% |
| Repeatability | Brix: ±0.1% |
| Dimensions & Weight | 3.2×3.4×16.8cm, 130g |

Options

- · 30% Sucrose (±0.03%) : RE-110030
- · 40% Sucrose (±0.04%): RE-110040
- · 50% Sucrose (±0.05%): RE-110050
- · 60% Sucrose (±0.05%): RE-110060
- · Daylight plate for small volume samplings : RE-2311-58M
- * Calibration Certificate: Contact an ATAGO representative for details,


Parts


· Daylight plate: RE-2315-60M

Design Registration No.000379326-0001, 000379326-0002 (EU) ZL 2005 3 0116403.4,5 (China), D111526 (Taiwan), D555, 549 (U.S.A.), 1255763, 1255764, 1255765, 1255766, 1255767 (Japan) Patent Granted in countries around the world.



@ Enlargement

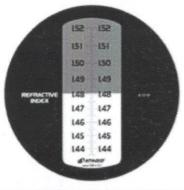

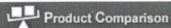







1.90MB - Download

MASTER-RI is refractive-index model (measurement range: RI 1.435 to 1.520) of the MASTE master-RI Series.

MASTER-RI can be used for measuring the refractive index of vegetable oil, petroleum, oil slicks, etc.

*Calibration with the Standard liquid LK (RE-99010)

MARTER PI

Specifications

| Model | MASTER-RI | |
|----------------------|----------------------------------|--|
| Cat.No. | 2612 | |
| Scale | Refractive index | |
| Measurement Range | Refractive index: 1.435 to 1.520 | |
| Minimum Scale | Refractive index : 0.001 | |
| Dimensions & Weight | 3.2 x 3.4 x 16.8cm, 130g | |

Options

- Standard liquid LK: RE-99010 - Daylight plate for small volume samplings : RE-2311-58M
- * Calibration Certificate: Contact an ATAGO representative for details.

Parts

- Daylight plate: RE-2315-60M

Honey refractometer (Automatic Temperature Compensation & Water Resistant) MASTER-HONEY

Design Registration No.000379326-0001, 000379326-0002 (EU) ZL 2005 3 0116403.4,5 (China), D111526 (Taiwan), D555, 549 (U.S.A.), 1255763, 1255764, 1255765, 1255766, 1255767 (Japan)

Patent Granted in countries around the world.

Tay PANA MAY EL

A TOP